


excellence
line

excellence line


VARIO magnetic holder

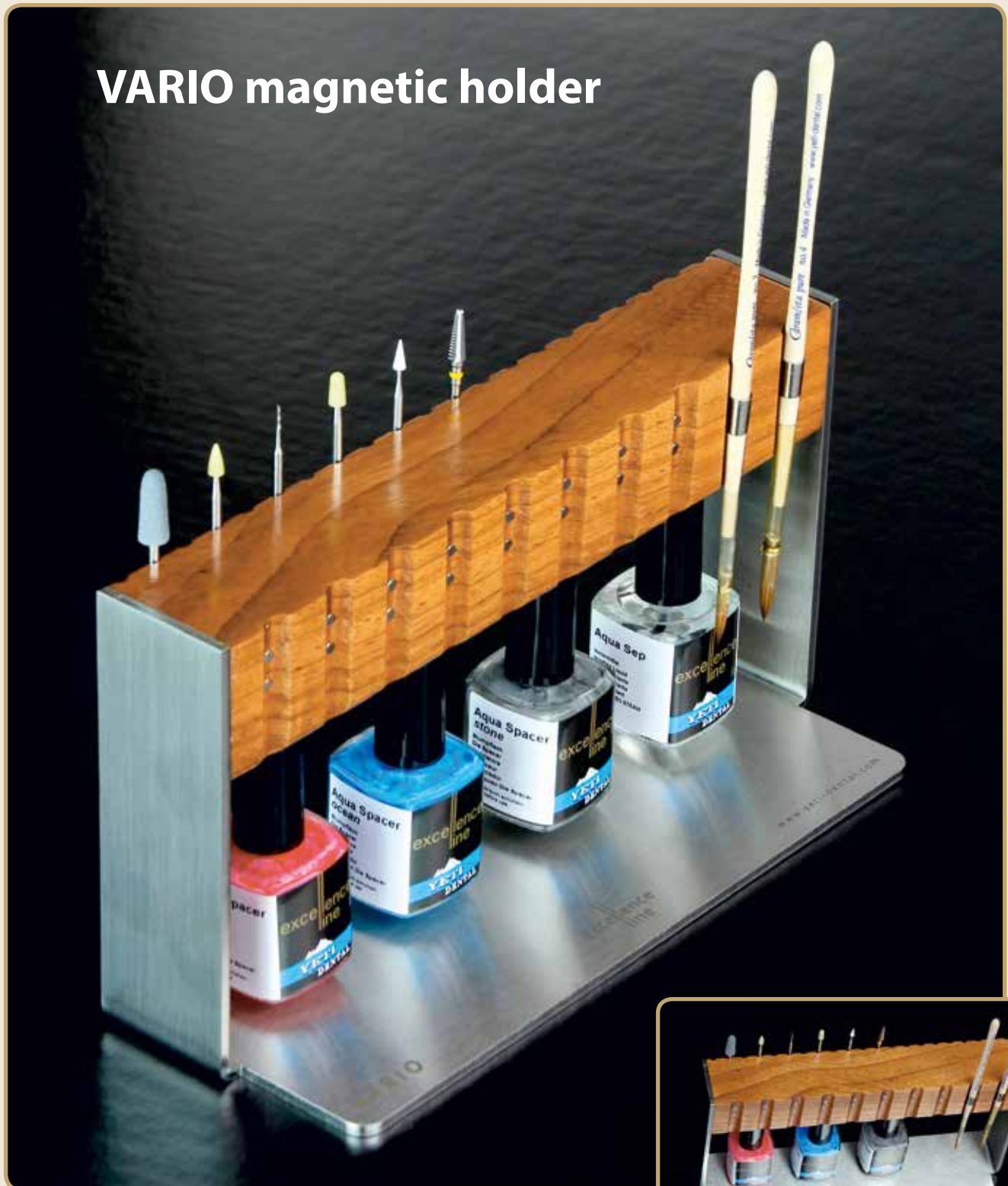

The Vario magnetic holder represents the centrepiece of our “excellence line”. Being made of elegant and finest materials, this all-rounder offers the possibility of magnetic fixation of our “excellence line” products. This is how Aqua Spacer, grinding tools, K2 Lube and Caramista brushes find their place on magnets immediately as if by magic. Your working place stays always tidy, all products are clearly arranged and immediately at hand. This saves place and time. Additionally, the high quality materials such as brush tips are protected due to hanging storage. Thus you can enjoy our excellence products for a long time.



Magnetic fixation of the “Caramista brushes” for optimal hanging storage
- Protection of the high-quality brush tips

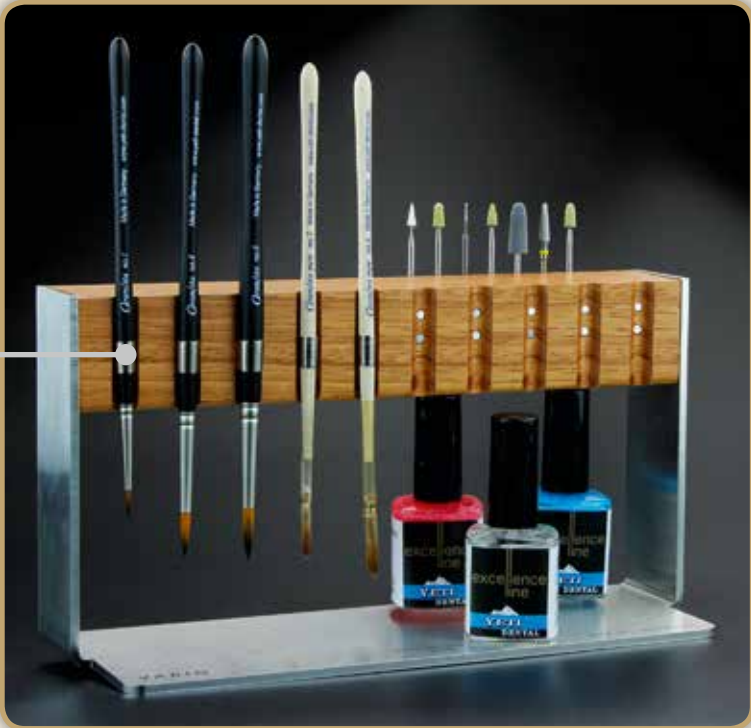

Lids of Spacer and K2 Lube fix to the Vario magnetic holder by magnets





Aqua Spacer, 15 my

Solvent-free spacer

Adheres well to the Vario Magnetic holder, always tidied up and close at hand!

- water-based
- can be easily diluted by distilled water
- solvent-free –no irritating vapours

Only use with Aqua Sep Separating liquid!

Content: 18 ml
Layer of 15 my



Aqua Spacer
coral



Aqua Spacer
ocean



Aqua Spacer
stone



Aqua Sep
Separating
liquid



K2 Lube

ceramic separating liquid

Suitable for the excellence series, it fixes magnetically to the Vario magnetic holder.

This ceramic separating liquid ensures reliable separation of ceramic from gypsum. It burns without residue and causes no colour changes of ceramic.

Content: 18ml



Caramista pure

Natural quill brush

The exclusive Kolinsky Hair is traditionally bound individually by hand and clamped in natural quill.

Natural hair brush from high-value Kolinsky sable hair:

- Perfect tension
- thick hair as water reservoir
- thin brush point

Brush handle:

- With metallic clip for fixation on the Vario magnetic holder

Comes in a natural linen brush pocket

- for protection of the brush hair
- with additional slots for other preferred instruments

„Handmade in Germany“

Content:

Set of Caramista pure no. 2 and no. 4 in handy brush pocket



All brushes of the series “Caramista” fix to the Vario magnetic holder by magnets. Up to 10 different brushes can always kept safely at hand.

Other brushes of series “Caramista” can be found in our main catalogue and at www.yeti-dental.com.

Caramista

Caramista N




excellence
line

„excellence line“
Zircon polishing set




Polishing set for the creation of natural texture of the tooth surface and for the production of zirconium frameworks

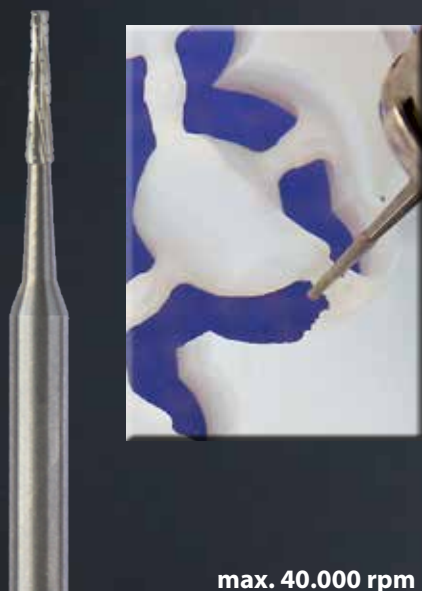

max. 40.000 rpm

HM Fissure bur
for separation of objects
from zirconia blank

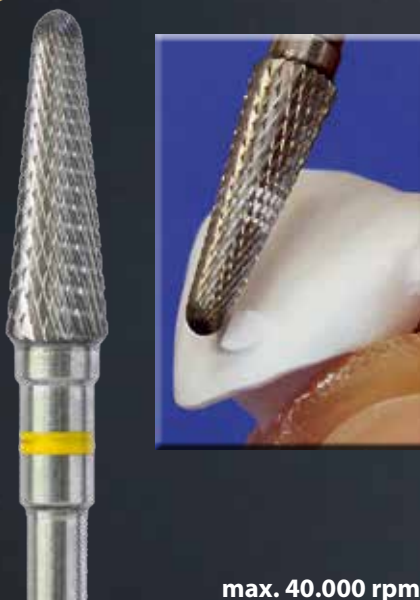

max. 40.000 rpm

HM cutter
in diamond form extra fine
for quick material removal
and very fine surface finishing
(suitable for layering ceramic)



max. 12.000 rpm

Smart Stone (2 forms)
with diamond grain for smooth
grinding of the contours
(suitable for sintered objects)



max. 20.000 rpm

Arkansas
for contouring of the surface
and and for finest details



max. 10.000 rpm


Abrasive paper cone
for smoothing of the surfaces
and imitation of the growth lines
(Perikymatien)



max. 10.000 rpm

Rubber Flame
for smoothing of the desired details,
for example abrasion surfaces

| Article | | Content | Article Number |
|---|---|---|------------------------------------|
|  | | Vario magnetic holder without instruments | 1 pc. |
|  | | Zircon polishing set 1 HM Fissure bur 1 HM cutter 1 Smart Stone cone 1 Smart Stone flame 1 Arkansas 1 Abrasive paper cone (mounted on mandrell) 1 Rubber Flame | 1 Set = 7 pcs. |
| Refill packages for zircon polishing set |  | HM Fissure bur | 5 pcs. |
| |  | HM cutter | 1 pc. |
| |  | Smart Stone cone | 1 pc. |
| |  | Smart Stone flame | 1 pc. |
| |  | Arkansas | 5 pcs. |
| |  | Abrasive paper cone (unmounted, without mandrell) | 25 pcs. |
| |  | Rubber Flame | 5 pcs. |
|  | | Aqua Spacer coral | 1 x 18 ml |
| | | Aqua Spacer ocean | 1 x 18 ml |
| | | Aqua Spacer stone | 1 x 18 ml |
| | | Aqua Sep Separating Liquid | 1 x 18 ml |
|  | | Caramista pure Brush Set | 2 brushes: 1 x no 2 1 x no 4 |
|  | | K2 Lube | 1 x 18 ml |

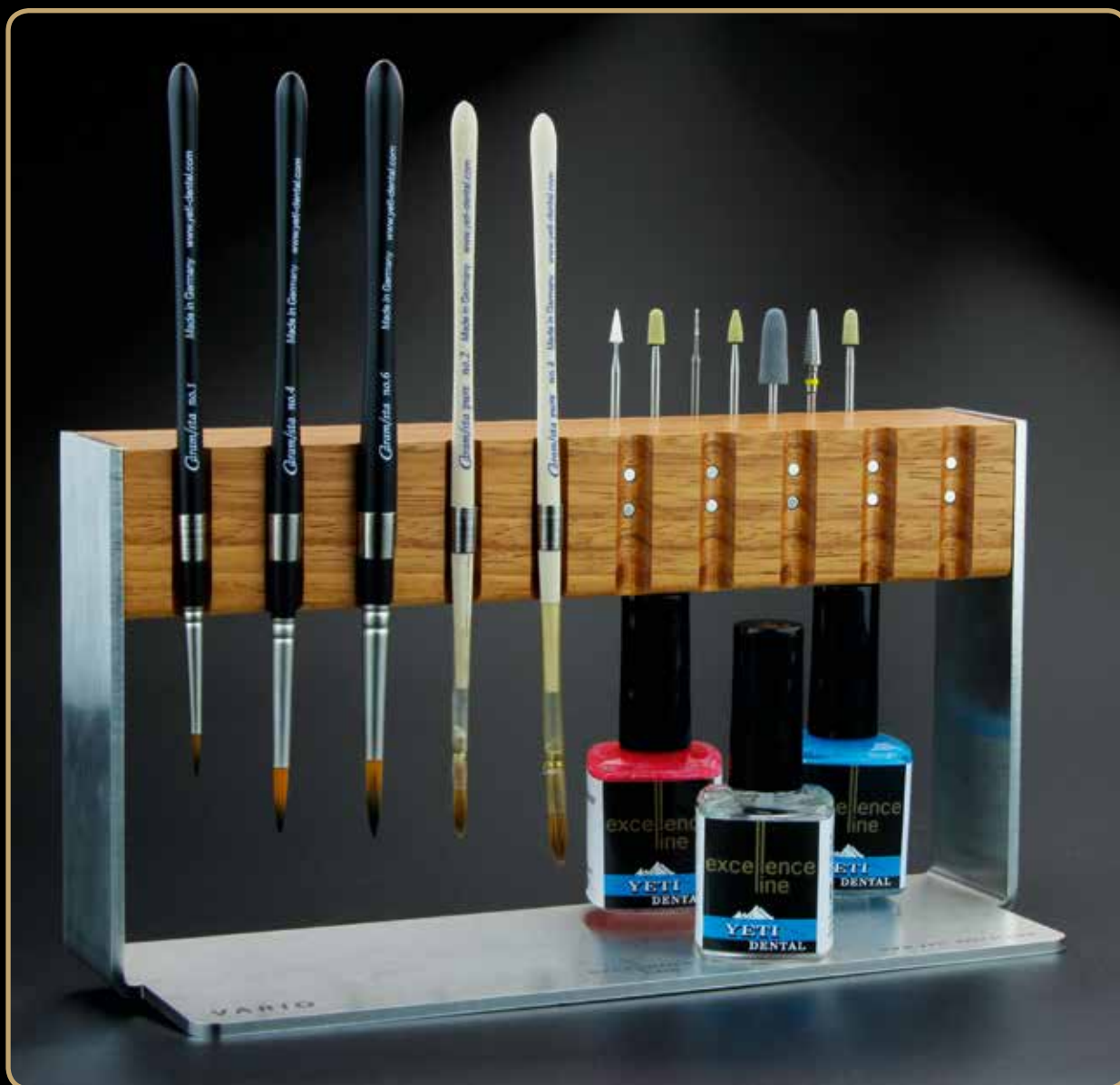

Print: BaurOffset Print e.K.

Konzept & Design: Anke Biskupski, Engen

03.2017 V1 by YETI Dental